POSITIVE POSTURE®

COREWAVE

OWNER'S MANUAL

Read this manual thoroughly before you start using your massage chair, paying close attention to the sections of this manual devoted to safety.

Follow the instructions in this manual for the safe use and maintenance of your massage chair. This manual should be kept on hand as a reference.

This product is designed for home use.

POSITIVE POSTURE
COREWAVE MASSAGE CHAIR

CONTENTS

02 Introduction

03 Cautions & Warnings

11 Main Features

13 Parts & Functions

15 Assembly

19 Moving the Chair

20 Sessions

22 Operations & Use

35 Troubleshooting

38 Warranty & After Sales Service

39 Specifications

Note: The pictures in this manual are only for reference. Please refer to the actual product

INTRODUCTION

Welcome to the world of Positive Posture and the transformative experience of your new CoreWave massage chair. More than just a luxury, this chair serves as your personal sanctuary for restoration and deep physical renewal.

Engineered with our ground-breaking CoreWave technology, this chair delivers a specialized combination of targeted vibration therapy and precise massage that revitalizes muscles in ways no other massage chair can. Combined with our advanced dual mechanisms, the CoreWave offers a sophisticated, synchronized massage experience that works multiple muscle groups simultaneously for unparalleled efficiency and relief. We have crafted this chair to be your dedicated partner in wellness, empowering you to prioritize your self-care and face each day with a renewed sense of vitality. We are confident that the CoreWave will become an indispensable part of your daily routine, and we sincerely hope you cherish the exceptional care and lasting benefits it provides.

CAUTIONS & WARNINGS

Indications, Symbols, & Descriptions

 **WARNING**	This action could result in serious injury or death
 **CAUTION**	This action could result in personal injury or damage to property
	This action is prohibited
	Further action is required
	Do not disassemble
	Unplug power cord from wall outlet

Important Safety Precautions

- Please follow all basic precautions while using this chair.
- Please read all instructions in this Owner's Manual before assembling and using the massage chair as it contains important information about assembly, and correct, safe procedures for using the chair.
- The above symbols are used to identify important safety instructions. Please familiarize yourself with these symbols and their meanings.
- Remember to keep this Owner's Manual in an easily accessible place.

General Safety



TO REDUCE THE RISK OF ELECTRIC SHOCK

Always power off and unplug the chair from the electrical outlet after use and prior to cleaning.



TO REDUCE RISK OF BURNS, FIRE, ELECTRIC SHOCK, OR INJURY TO PERSONS

- Only plug the chair into a properly grounded outlet. See Grounding Instructions on page 5.
- The chair should never be left unattended when plugged in. Turn off and unplug from outlet when not in use, befor cleaning, moving, or making physical adjustments to the chair.
- Do not operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- Close supervision is necessary when the chair is used by, on, or near children or individuals with disabilities.
- Use the chair only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate the chair if the cord or plug is damaged, if it is not working properly, if it has been dropped or damaged, or if it has been exposed to water. If this happens, contact the service center for examination and repair.
- Do not carry the chair by the supply cord or use the cord as a handle.
- Keep the chair and its components away from heated surfaces.
- Never drop or insert any object into any opening.
- Do not use outdoors.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- To disconnect, turn all controls to the off position, then remove plug from the outlet.
- Use heated surfaces carefully. May cause serious burns. Do not use over areas with sensitive skin or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
- Keep children away from extended foot support (or other similar parts).
- Do not use the chair in an excessively humid or dusty environment as this may result in malfunction or electric shock.

CAUTIONS & WARNINGS

Grounding Instructions

This product must be grounded. If it should malfunction or break down, grounding provides a path of least resistance for electric shock. This product is equipped with an equipment-grounding conductor and a grounded plug. The plug must be plugged into an appropriate outlet that is properly installed and grounded in accordance with all local codes and ordinances.

DANGER: Improper connection of the equipment-grounding conductor can result in a risk of electric shock. Check with a qualified electrician or serviceman if you are in doubt as to whether the product is properly grounded. Do not modify the plug provided with the product—if it will not fit the outlet, have a proper outlet installed by a qualified electrician.

This product is for use only on a nominal 120V circuit and has a grounded plug that looks like the plug illustrated in Figure A below. Make sure that the product is connected to an outlet having the same configuration as the plug. An adapter should never be used with this product, such as shown in Figure B below.

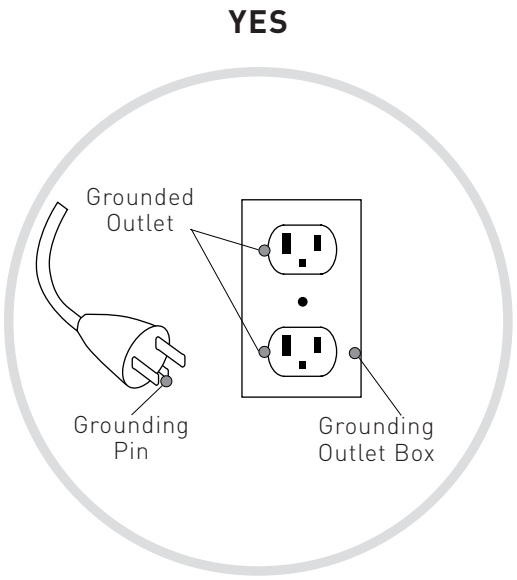


Figure A

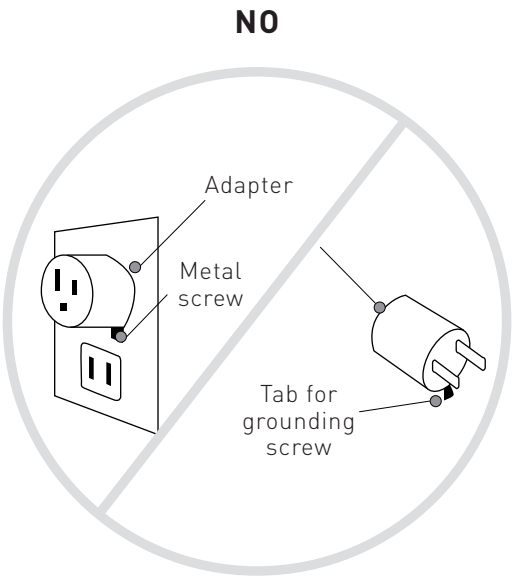


Figure B

Safety During Installation

WARNING

ELECTRICAL GROUNDING

- This product uses a normal Alternating Current (AC) circuit and has a ground pin as displayed in *Figure A* on page 5. Always make sure that the chair is connected to an outlet that has the same configuration as the plug, and is fully inserted into the outlet. If the chair should malfunction or break down, the grounding will provide a path of least resistance for electrical current to help reduce the risk of electric shock.
- This chair must be plugged into an appropriate outlet that has been properly installed and grounded in accordance with all local codes and ordinances. We strongly recommend using a high-quality surge protector with your chair. A surge protector offering at least 2,000 joules of surge suppression is required for warranty service.
- Do not use with any type of transformer device.
- Do not use the chair with a power supply other than those shown in *Figure A* on page 5.
- Make sure the main switch on the back base of the chair is in the OFF position before inserting the plug into the outlet.
- Any modification of the plug provided is prohibited and can cause serious injury. If the plug does not fit into an outlet, contact a qualified electrician for proper installation.
- Failure to connect the grounding conductor of the chair can result in the risk of electrical shock. If you feel the product is not properly grounded, please cwontact a qualified electrician.
- Please unplug this chair from the electrical outlet after use and before any cleaning or maintenance. To do so, ensure the main power switch is in the OFF position, then unplug the chair from the outlet.

CAUTION

ELECTRICAL GROUNDING

- This product is intended for use described in this Owner's Manual. Do not use any attachments that are not recommended by the manufacturer.
- Do not use the chair if the cord or plug is damaged, if it is not working properly, if it has been dropped or damaged, or if it comes into contact with water. In case of any of these scenarios, please contact Positive Posture Service immediately.
- Make sure all air openings are free of lint, hair, or other debris. Do not operate the chair with the air openings blocked.
- Never use the chair where aerosol spray products are used or oxygen is being administered.
- Excessive heating of this product can occur and cause fire, electric shock, or serious injury. Do not operate under a blanket, pillow, or heated surface.
- Do not try to move or carry the chair by pulling or holding onto the power cord.
- Keep power cord away from all heated surfaces.

CAUTIONS & WARNINGS

Safety During Installation



TO AVOID ACCIDENT AND INJURY

- Install the chair on a flat and level surface.
- Allow at least 6" (15 cm) clearance from the chair's back to a wall to allow sufficient space for reclining up and down.
- Do not drop or slam the footrest.
- This chair is not recommended for children under the age of 14. Please keep children away while in use or storage.

TO AVOID DAMAGE

- Do not use outdoors.
- Ensure the power cord is not twisted or pinched by any part of the chair or surrouding items/furniture.
- Stoves, heating products, direct sunlight, or other high temperature areas can cause fading, discoloration, or hardening of the upholstery. Please keep your chair away from these types of exposures.

Safety Before Use



PEOPLE NOT SUITABLE TO USE THE CHAIR

- People with circulatory problems (such as varicose veins or thrombosis), children under the age of 14, and those advised not to use a massage chair by their doctor should not use this chair.
- If you have any of the following health conditions, please consult your doctor for medical advice before using this chair: serious heart problems, serious diabetes, serious osteoporosis, serious skin conditions, malignant tumor(s), sensory impairment, pregnancy (or suspected pregnancy) or those who have just given birth, acute gastrointestinal complaints (gastritis, hepatitis, enteritis), joint dysfunction such as rheumatism, hammer toe and gout, thecitis (or suspected thecitis), high fever, spinal abnormalities due to past injury or ailment, curvature of the spine (scoliosis), acute neck sprain (whip lash), or any other medical condition that may contraindicate usage of this product.
- Seek medical advice before using this chair if you have a pacemaker or other electronic medical device, if you have been restricted to bed rest, or if you are planning to use this chair to treat a specific ailment, recover from surgery, or as physical therapy.
- Elderly people and those with weak bones should seek medical advice before using this chair even if they do not have a specific disease or disorder mentioned above.
- Do not use this chair if impaired by alcohol or other substances.
- The chair is not intended for use by persons with reduced physical, sensory or mental capabilities, or who lack experience and knowledge of proper chair operation, unless they have been given appropriate instruction concerning use of the chair or are under the direct supervision of a person responsible for their safety.

Safety Before Use, continued...



TO AVOID ACCIDENT, INJURY, AND DAMAGE

- Children should be supervised to ensure that they do not play with the appliance.
- The appliance has a heated surface. Persons sensitive to heat must be careful when using the appliance.
- Do not use this chair for any other reason than described in this manual.
- Do not use this chair with any other therapeutic device or accessory besides those included.
- This chair is intended to be used with all the fabric attached. Failure to do so can result in serious injury.
- Check for any holes or tearing in the active massage area of the chair before each use. If a tear is noticed, stop using immediately and contact Positive Posture Service.
- Check the position of the roller mechanism prior to sitting down. If the rollers are not in the storage position, turn the power on and check that the chair is operating properly. For the proper storage position, see page 14.
- Do not move the chair while someone is in it.
- Do not stand on this chair.
- Do not sit on any part of the chair other than the seat.
- Do not sit, stand, or lean on the footrest unless you are sitting in or getting up from the chair.

Safety During Use



TO AVOID ACCIDENT AND INJURY

- Do not wrap the power or remote cord around your body or place the chair on top of the cord.
- Do not scratch, tear, treat, twist, stretch, or bend the power or remote cord.
- Spilling water or other liquids on this chair may result in electric shock, short circuit, or combustion.
- Keep all children and pets off the chair and away from the surrounding moving area.
- Do not allow children or individuals with disabilities near the chair without close supervision. Keep all persons away from the reclining backrest and extending footrest to avoid injury.
- Do not use the chair while more than one person is sitting in it, or when holding a baby or small child.
- Do not use the chair if any part of your body, including your hands, is wet.
- Do not use the chair unclothed.
- Do not use this chair for any purpose other than its intended uses outlined in this manual.
- Do not use the chair to massage the head, chest, stomach, and/or joints including knees and elbows.

CAUTIONS & WARNINGS

Safety During Use, continued...



TO AVOID ACCIDENT AND INJURY

- Do not place your hands or feet into the roller tracks, or into any moving parts, such as the reclining mechanism.
- Remove all hard objects, such hair accessories or jewelry, before using.
- Do not use the leg rest for anything other than legs or feet.
- Do not operate this massage chair without someone in it.
- Do not put your hands or head between or beneath the leg rest.
- Do not sleep in the chair.
- Never unplug the power cord during operation.
- Upon initial use of this chair, start off with a lighter massage. Starting off with a deep massage may result in injury. Elderly and people with weak bones should be careful while selecting massage intensity.
- If you notice any abnormal condition or in case of an emergency, push the Pause button on the remote to stop operation immediately.
- Immediately discontinue use of the chair and seek medical advice if you notice any acute pain or unusual physical discomfort.
- If the massage feels too strong, reduce the stimulation according to this manual. If it still feels too strong, discontinue use immediately.



TO AVOID ACCIDENT AND INJURY

- Do not use this chair if there is a power outage or the possibility of one. In this case, stop using the chair immediately, turn the power switch off, and unplug the cord from the wall outlet.
- Do not use this chair if there is a threat of lightning. If so, stop operation immediately, turn the main power switch off, and unplug the chair from the wall outlet.
- If a specific operation does not start or you find any incidence of abnormality, stop using the chair, turn off the main power switch, and unplug the cord from the wall outlet. Call Positive Posture immediately.

Safety During Use, continued...



TO AVOID DAMAGE

- Do not use this chair for any longer than 30 minutes each day. Avoid using a specific focused area massage for more than 5 minutes at a time. Long and continuous massage in the same position may cause discomfort or bodily harm.
- Do not attempt to force your leg, foot, arms, or hands back into the leg rest/arm rest if it comes out. Undue force may result in damage to the chair.

OTHER SAFETY PRECAUTIONS

- Do not cause any major impacts to the chair.
- Do not sit on any part of the chair other than the seat.

Safety After Use and When Not in Use



TO AVOID ACCIDENT, INJURY, AND DAMAGE

- Do not plug or unplug the chair with wet hands.
- Do not use any other power cord except the one supplied.
- Grasp the power plug when unplugging the cord from an outlet. Do not use the power cord to pull.
- Turn the main power switch off and unplug the power cord from the wall outlet after each use.
- Always recline the seat back to the upright position after using.
- Wipe down the chair after each use.
- Always store the chair in a low humidity environment.
- Do not place any objects on the chair.
- In case of malfunction, do not attempt to repair. Contact Positive Posture Service.
- Do not modify any part of the chair.
- Be sure to include this Owner's Manual when transferring this chair to another party.
- Do not dispose of this product or its accessories without consulting with your local authorities beforehand.

MAIN FEATURES

PROPRIETARY COREWAVE™ ULTRA HERTZ VIBRATION

Proprietary CoreWave™ Ultra Hertz vibration is a breakthrough in deep tissue massage therapy. The 4D upper massage rollers vibrate at a high frequency of 4,200/minute to melt away stubborn knots with pinpoint precision. This advanced technology penetrates deeper into the muscle tissue than conventional rollers, providing unparalleled, profound relief.

DUAL MASSAGE MECHANISMS

Experience total-body harmony with our synchronized Dual Massage Technology. Two independent mechanisms work simultaneously to melt away tension from top to bottom. The advanced upper 4D rollers (with proprietary CoreWave™ Ultra Hertz Vibration) masterfully target your neck, upper back, and lumbar, while a dedicated lower system delivers deep, soothing relief specifically to your sacrum and gluteal muscles for a truly comprehensive massage.

SL-TRACK

CoreWave's SL-Track guides the Dual Massage Mechanisms from neck to glutes for a comprehensive full-body massage.

INTELLIGENT BODY SCAN

Advanced body scanning technology intelligently maps your unique shape and size at the beginning of each session, locating acupressure points and tailoring each massage session to your body.

18 AUTOMATIC PROGRAMS

Discover the perfect massage with CoreWave's 18 expertly designed automatic programs. From full-body relaxation to targeted relief, each program offers a unique combination of massage techniques to suit your needs.

IMMERSIVE FULL BODY AIR MASSAGE

Strategically located air cells provide customizable levels of pressure to mobilize your body, improve circulation, and deliver targeted relief to your shoulders, arms, hands, hips, legs, and feet.

DUAL ZONE HEAT THERAPY

CoreWave's Dual Zone Heat Therapy offers three heat levels to melt away tension, fatigue, and stiffness as it prepares muscles for massage therapy by applying targeted heat to the calves and lower back.

JAPANESE-STYLE FULL EXTENSION STRETCH

Fully mobilize and decompress your lower back and lumbar as CoreWave elevates your core while stabilizing your shoulders and legs.

LOWER LEG TRACTION

CoreWave's footrest extends and retracts to gently decompress your ankle joints, relieve built-up tension, and restore flexibility after a long day on your feet.

TRIPLE ACTION REFLEXOLOGY FOOT ROLLERS

Soothe tired feet and stimulate wellness-enhancing pressure points with the touch of three expertly-crafted mechanisms on each foot that individually address the heels, arches, and balls of your feet.

AUTOMATIC LEG EXTENSION

CoreWave uses intelligent sensors to instantly measure your leg length. The foot rest automatically extends or retracts to give you a custom fit, maximizing your overall comfort and ensuring perfect alignment with the calf massagers.

THERAKNEAD™ CALF MASSAGE

Soothe your calves with a masterful combination of mechanical massage, gentle kneading and therapeutic compression. Air cells compress and rotate in a graceful circular motion while massage rollers apply precise pressure, melting away tension, revitalizing tired muscles and encouraging healthy circulation—leaving your legs refreshed, renewed, and ready to carry you forward.

TRUE ZERO GRAVITY® RECLINE

Experience weightless relaxation in True Zero Gravity® recline. CoreWave offers three preset positions that elevate your legs and evenly distribute your weight, minimizing pressure on your spine and promoting deep relaxation. Float away as you recline into perfect equilibrium.

MAIN FEATURES

MEMORY SETTINGS
Save up to six of your auto-program massage settings, including recline position, intensity, heat level, and air massage compression level, so your perfect massage is available with a simple touch of a button.

TOUCHSCREEN CONTROLLER
Effortlessly control your massage experience with the intuitive touchscreen remote. Easily navigate through massage programs, adjust intensity levels, and target specific areas with a simple touch.

ADJUSTABLE TIMER
Whether you need to get started a bit early to seize the day or stay for a few extra minutes of therapy, easily adjust the duration of your massage.

BLUETOOTH® SPEAKERS
Relax to your own soundtrack with CoreWave's built-in Bluetooth® sound system.

ADJUSTABLE HEADREST
Designed for personalized comfort, the adjustable headrest pillow perfectly cradles your head, so you can find the perfect position and melt into pure bliss.


REMOVABLE BACKPAD
Further customize the intensity of your massage with the removable back pad. Keep the pad on to enjoy a gentler, more soothing massage, or remove the pad for a more intense, therapeutic massage.

INTELLIGENT VOICE COMMAND
Activate your favorite massage programs, adjust your recline angle, and more, without lifting a finger.


AMBIENT LIGHT
Enhance your relaxation with warm ambient light that gently fills the room.

SPACE SAVING DESIGN
CoreWave glides forward as it reclines, requiring only 6" of distance from the wall when upright. This clever feature allows you to enjoy a full body massage without sacrificing valuable living space.

PARTS & FUNCTIONS

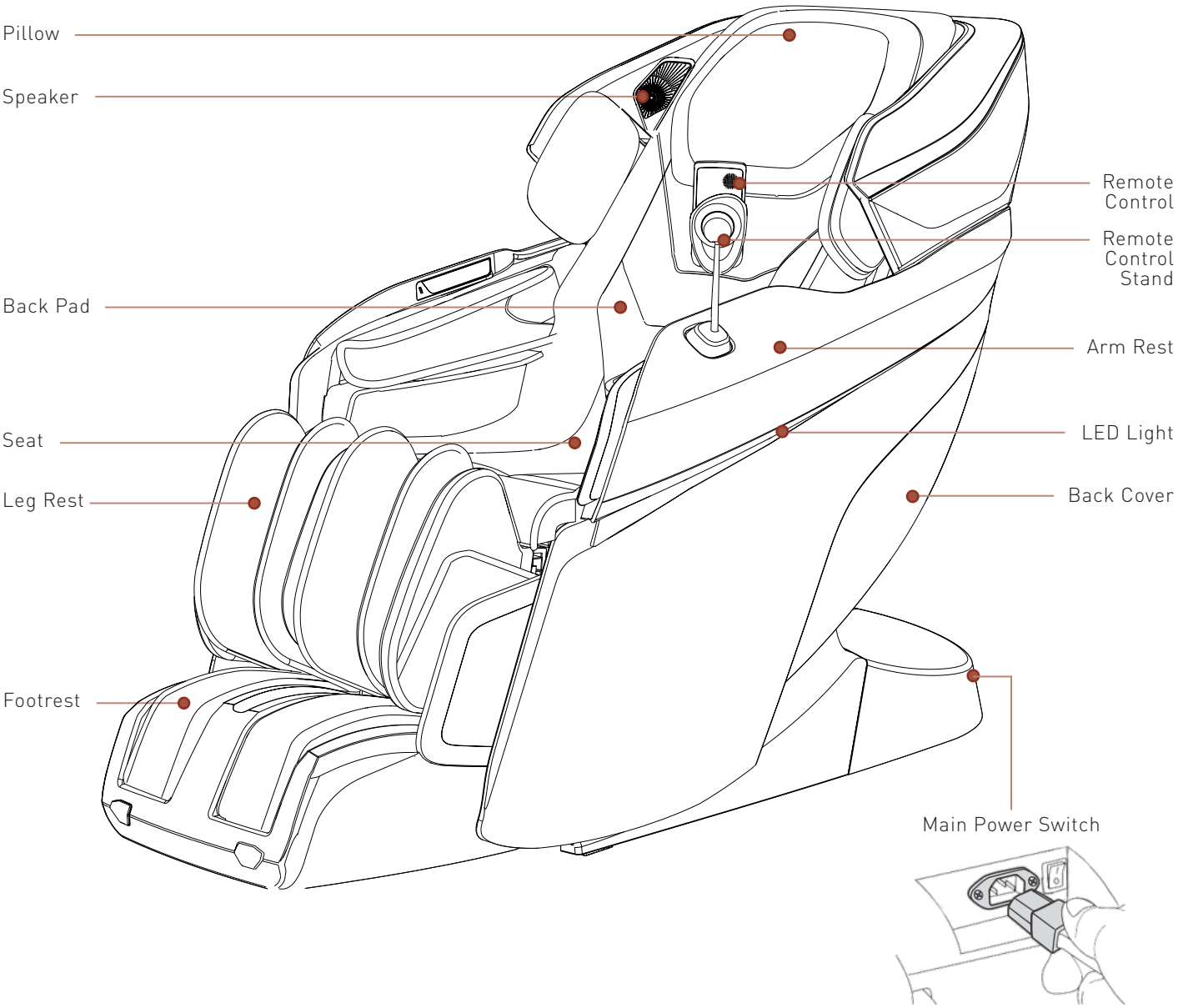

PARTS & FUNCTIONS

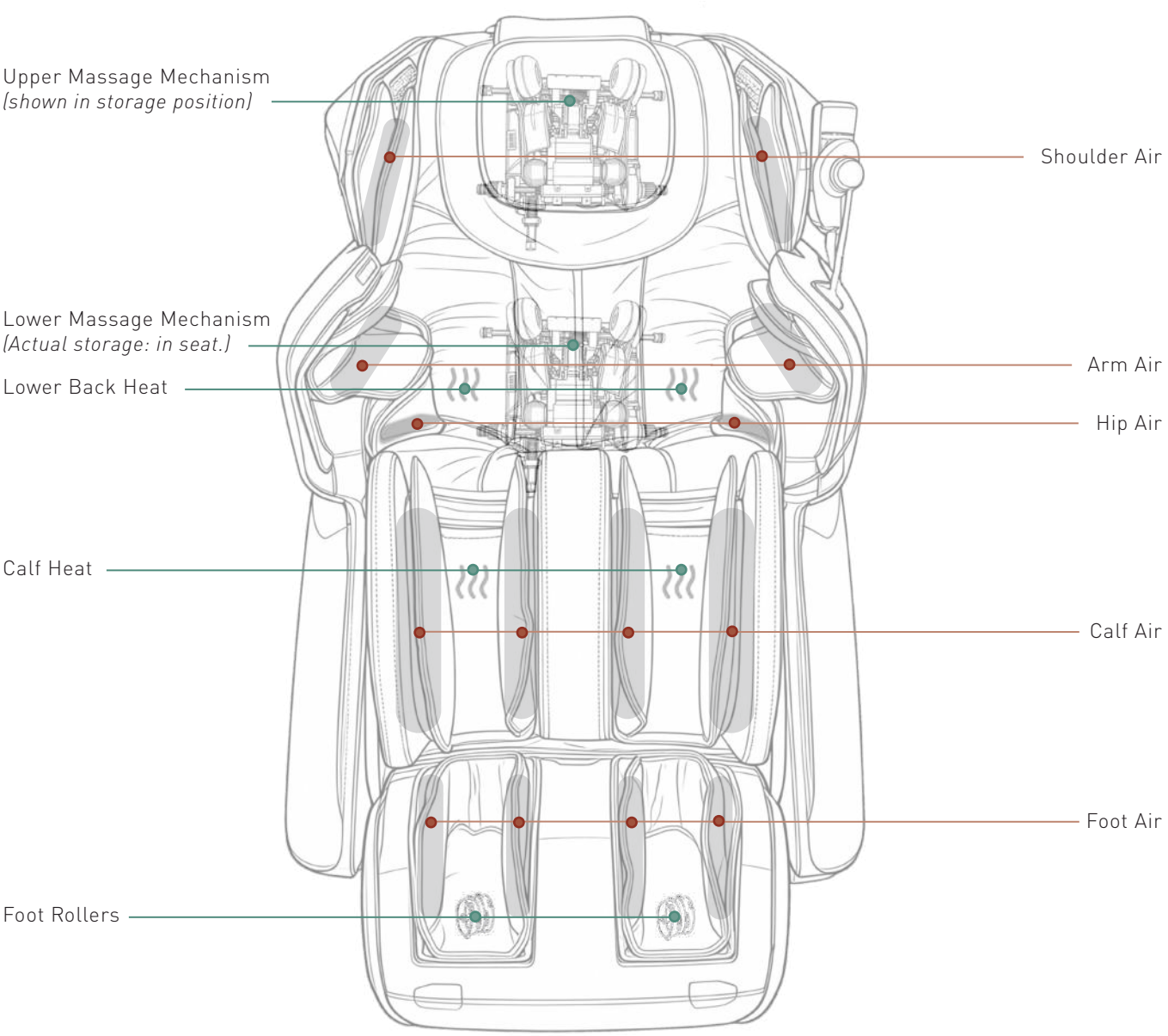

ASSEMBLY

Safety During Installation

CAUTION

TO AVOID DAMAGE

- Do not install the chair in an area of high humidity, such as a bathroom or sauna.
- Do not use the chair in a room with a temperature of 104-degrees or higher.
- Do not expose the chair to heaters, stoves, or direct sunlight.
- Keep fire sources, such as cigarettes and ashtrays, away from the chair.
- Ensure chair is installed on a flat and level surface to avoid tipping over.
- As the chair is heavy, please exercise extra care not to damage the floor when positioning the chair.

Assembly Preparation


- Ensure there is an allowance of at least 6" (15 cm) behind, and 15.75" (40 cm) in front of the chair to allow the backrest to recline and the leg rest to extend without obstruction during operation.
- Maintain a minimum distance of 3' (92 cm) from AV devices, such as a TV or radio, to prevent signal interference.
- It is recommended that a carpet be placed under the chair to prevent scratching the floor and to reduce noise during operation.
- As the chair is heavy, please exercise extra caution.




Base




Leg Rest

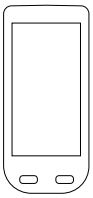


Remote

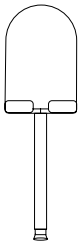

Remote Stand

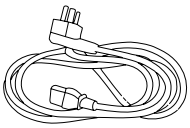

Power Cord

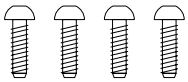

4 Screws



Leg Rest Mounting Bracket (2)



Fuse (*Non-Assembly Spare Part)*

Assembling the Chair

CAUTION


- To prevent disconnection, be sure to double check all connections by gently pulling on the electrical connectors and air hoses to ensure they are securely attached.
- If not connected properly, the air cells will not inflate as designed.

STEP 1

Unbox chair and place in designated area. Carefully remove all items (including leg rest from the seat).

STEP 2

Line up the remote control stand with the slot on the left armrest and twist in until secure. Then, place the remote control on the remote control stand.

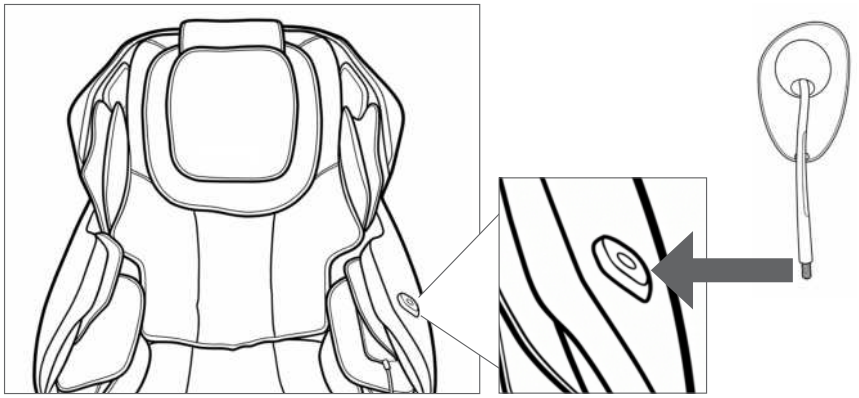

ASSEMBLY

STEP 3

- Line up the leg rest with the front base of the chair, then attach the electrical line and the air line from the leg rest to the chair's base (*Figure 3A*).
- Slide the calf unit mounting brackets onto the left and right side posts found on the top of the leg rest [*Figure 3B*].

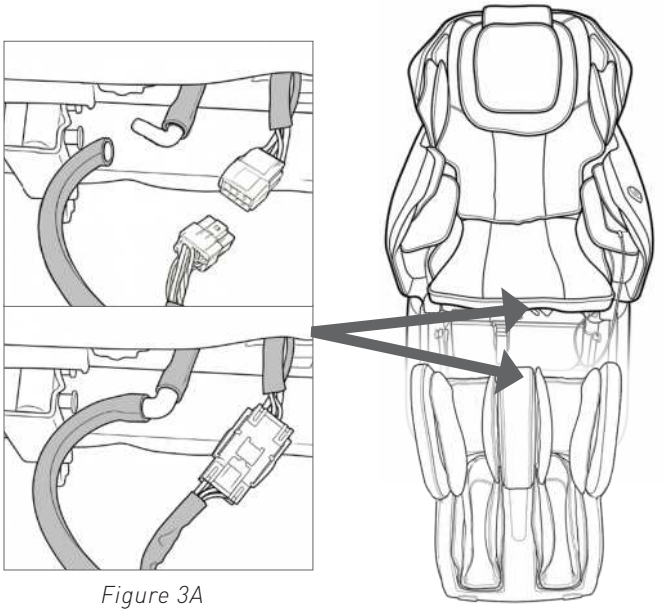


Figure 3A

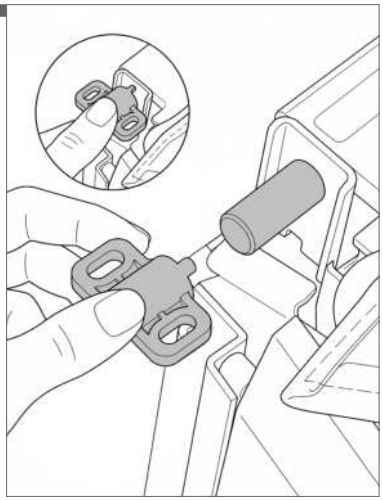


Figure 3B

STEP 4

- Secure the leg rest by carefully lifting and placing the calf unit mounting brackets on the top of the leg rest onto these brackets *(Figure 4A)*. Be gentle to avoid damaging the electrical connection.
- Install two screws into each bracket to secure *(Figure 4B)*.

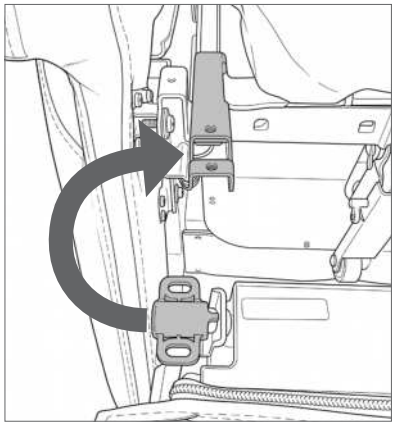


Figure 4A


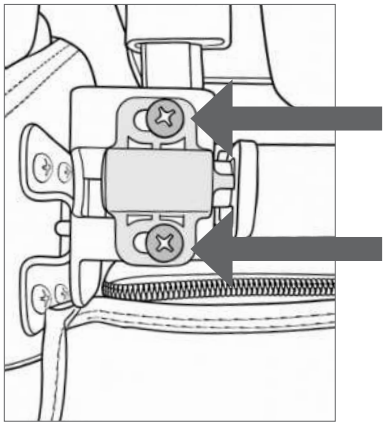


Figure 4B


STEP 5

Carefully place the back/seat pad onto the chair's seat. Secure by connecting the zipper from the base of the chair to the seat pad on both the left and right side.

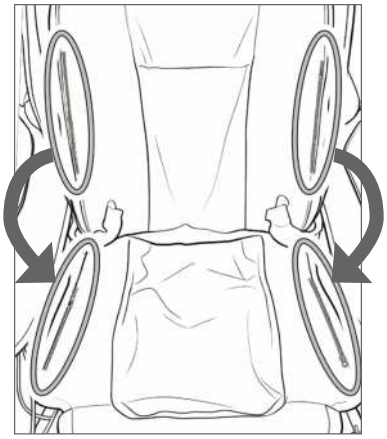

STEP 6

Secure the zipper from the front of the seat cushion to the top of the leg rest.

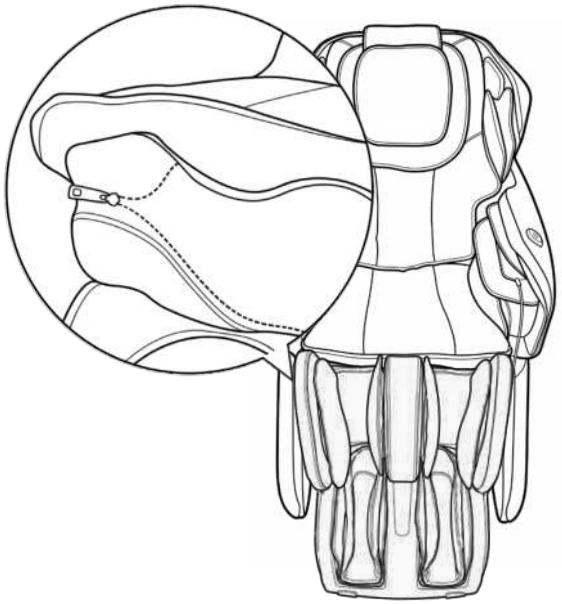

STEP 7

Turn the Main Power Switch on the back base of the chair to the ON position, then power on the remote control and begin using the chair.

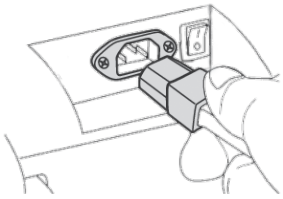


Figure 7A

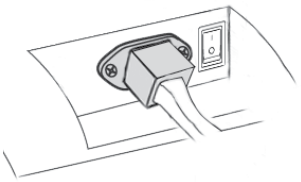


Figure 7B


MOVING THE CHAIR

CAUTION

TO AVOID DAMAGE

- Never move the chair while it is operating or when someone is in it.
- As the chair is heavy, please exercise extra care not to abruptly drop it. This may cause damage to the floor as well as the internal components of the chair.

Before moving the chair, turn off the main power and disconnect the power plug from the wall outlet. Bind the cables and put them in the seat to avoid damage when moving the chair.

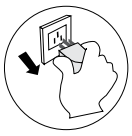


Option 1

When carrying the chair, the leg rest should be removed prior to lifting. Have one person carefully lift from the front of the chair while a second person lifts using the handle on the back of the chair. Carry to desired location.



Option 2

To move the chair without lifting, have one person press down firmly on the headrest while a second person lifts from under the bottom of the armrest to shift the weight onto the rear castors. Once the rear castors are engaged, carefully push the chair to your desired location, then gently lower it to the floor.

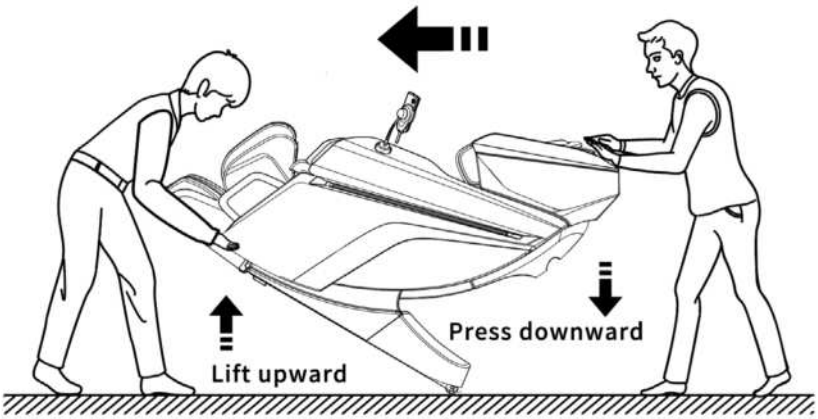

Programmed Sessions

Choose from eighteen programmed and focused massage sessions. Each session lasts between 10 and 20 minutes and can be selected by navigating through the easy-to-use menus on the remote display.

FULL BODY	Experience a complete physical and mental reset with the Full Body massage. This comprehensive experience wraps your entire being in a blanket of comfort, methodically soothing you from head to toe. It is the perfect choice when you want to escape the noise of the world and emerge feeling wholly balanced and completely renewed.
REFRESH	There's no better way to end a long day! Let the crisp, clean energy of the Refresh session provide the mental reset you need. This session is designed to clear the mental fog and lift a heavy mood. You'll step out of the chair feeling instantly lighter, brighter, and full of vibrant energy.
DEEP SHIATSU	Honor your body's need for profound, deliberate relief with Deep Shiatsu. This very active session mimics the pressure and movements of TWO professional therapists working on you at the same time. It will help you work away the deep-seated tension you've been carrying, leaving you with an empowering sense of release and structural freedom. Make it your go-to when you need to unlock stubborn physical stress.
STRESS RELIEF	Melt away the demands of modern life with the sanctuary of the Stress Relief massage. Specifically tailored for those days when your mind is racing and your shoulders are up by your ears, this program acts as an emotional exhale. It gently coaxes your nervous system into a state of quiet safety, leaving you wrapped in a cocoon of calm.
STRETCH	Open up your body, extend your spine, and expand your mental boundaries with the liberating Stretch program. This session beautifully elongates your posture and opens your chest, mimicking the elegant, fluid release that come with carefully planned back extension movements. It is the ultimate antidote to a sedentary day, leaving you feeling delightfully limber, taller, and filled with a sense of promise as you look up and toward the horizon.
NECK & SHOULDER RELIEF	Shed the weight of the world with the focused Neck & Shoulder Relief massage. This session cradles and comforts the exact areas where digital fatigue and daily responsibilities accumulate most. You will feel movements targeting tightness, leaving your upper body feeling unburdened, looser, and better aligned.
NIGHT	Clear all thoughts and slide into a state of sweet slumber with the peaceful Night program. This deeply comforting session acts as a gentle bridge between a hectic day and a restful evening, slowing down your internal rhythm. It coaxes your body to let go of the day's final worries, leaving you profoundly relaxed and perfectly prepared to drift into a deep meditation or deep sleep.
AWAKE	Ignite your morning or conquer the afternoon slump with the motivating Awake program. This dynamic session acts as a catalyst for your day, stimulating your senses and stirring your inner drive. You will leave the chair feeling fully alert, sharp, and brimming with the enthusiastic momentum needed to take on whatever lies ahead.
WARM UP	Prepare your body and mind for peak performance or a busy day with the inviting Warm Up session. Designed to gently stoke your internal fire, this program coaxes your muscles out of hibernation and invites fluid mobility back into your joints. It leaves you feeling supple, primed, and enthusiastically ready to transition into physical activity.

SESSIONS

COOL DOWN	Transition gracefully from high gear to a state of peaceful recovery with the Cool Down massage. Perfect after an intense workout or a high-pressure presentation, this program gently lowers your tempo and rewards your body for its hard work. It leaves you with a soothing sense of closure, leaving your muscles feeling relaxed and your mind quietly content.
MEDITATE	Journey inward to a place of profound stillness with the slow and exquisite Meditate program. This session provides a deeply comforting, rhythmic backdrop that allows you to detach from external chaos and anchor yourself in the present moment. It leaves you in a state of serene mindfulness, wrapped in a blanket of absolute inner peace and clarity.
ROCK & RELAX	Surrender to the comforting, nostalgic sway of the Rock & Relax program. By mimicking a gentle, rhythmic cradling motion, this session bypasses a busy mind and speaks directly to your body's need for safety and rest. It is the ultimate escape for pure comfort, leaving you feeling deeply grounded, peaceful, and blissfully content.
LYMPHATIC FLOW	Purify your system and invite a feeling of inner lightness with the Lymphatic Flow massage. This specialized session focuses on a sweeping sense of movement, helping you shed the feeling of physical heaviness and fluid retention. Emerge feeling cleansed, balanced, and unburdened from the inside out.
CRADLE	Find ultimate sanctuary and emotional comfort in the tender embrace of the Cradle program. This unique massage envelopes you in a nurturing hold, providing a safe haven to completely let go of your defenses. It is perfect for those times you feel emotionally drained, leaving you feeling profoundly protected, comforted, and deeply cared for.
SLEEP	Escape into the cozy stillness of the Sleep program. Inspired by quiet, mesmerizing beauty, this session encourages your mind and body to slow down, disconnect, turn inward, and find warmth. It leaves you in a state of deeply peaceful hibernation, especially on those nights when you can't find the stillness necessary to get to sleep (or stay asleep). Can't sleep? Roll out of bed and choose Sleep.
CIRCULATION	Awaken your body's natural vitality and life force with the invigorating Circulation program. This session focuses on flow throughout your entire system, sweeping away fatigue from your extremities. You will step out of the chair feeling tingly, revitalized, and filled with a radiant, healthy glow.
LEGS & FEET	Ground yourself and restore your foundation with the targeted Legs & Feet massage. Perfect after a grueling day on your feet or a long journey, this session breathes new life into heavy, aching lower limbs. You'll experience an immediate sensation of weightlessness and relief, leaving you feeling ready to walk on air.
WAIST & GLUTES	Unlock your body's core power and release hidden tension with the Waist & Glutes program. This deep, stabilizing session focuses on the pivotal center of your body, melting away the tight grip of sitting or stress. You will leave the chair with a profound sense of decompression, lower-body freedom, and effortless physical balance.

Manual Sessions

This chair has 6 manual massage methods from which you can choose using the upper mechanism: Knead, Tap, Knock, Knead & Tap, Shiatsu, and Roll.

This chair has 5 manual methods from which you can choose using the lower mechanism: Knead, Tap, Knead & Tap, Knock, and Roll.

Manual Air Sessions

This chair has 3 manual air sessions from which you can choose: Full Body, Upper Body, and Lower Body.



TO AVOID INJURY

- Always check your surroundings before beginning a massage to ensure that no people, pets, or objects will be hit or squeezed by the reclining back or leg rest.



TO AVOID DAMAGE AND INJURY

- This massage chair is meant to be used for indoor household use only. The chair is not designed for commercial purposes.
- If you use the chair in a cold room, do not increase the temperature abruptly. Instead, increase the temperature gradually to a normal level.
- If the chair is stored in a cold place and is moved into a warm environment, wait one hour before use. Otherwise, performance may suffer, or the chair may malfunction due to condensation on the chair.
- Do not sit on the chair with the leg rest raised. Sudden heavy weight on the leg rest could cause damage to the chair, or it may result in accident or injury.
- Ensure the backrest returns to the upright position before sitting.
- Do not use the arm or shoulder unit for any area other than your arms and shoulders.
- If your arm comes out of the arm or shoulder unit or your leg/foot come out of the leg rest during a massage, do not try to force it back in.
- Do not massage your elbows.
- Do not put anything other than legs and feet in or on the leg rest.



TO AVOID DAMAGE

- Do not sit, step, or place heavy objects on the remote control.
- Do not drop the remote control.
- It is recommended that the remote control be placed on the remote stand when not in use.
- Do not exert excessive force on the remote control.
- Do not pull the remote control cord forcefully.
- Do not operate the remote control using sharp edged or hard objects.

OPERATION & USE

Positioning Yourself in the Chair

- Sit back fully in the chair so your entire back is against the backrest and your head is leaning into the pillow.
- Adjust the pillow so it sits just above your shoulders and at the base of your neck. If the pillow is positioned too high or too low, the intelligent body scan will not work correctly. *(Figure A)*
- Place you arms and elbows between the arm unit air cells with your palms facing down. *(Figure B)*
- Position your feet in the leg rest so your heels and rear calves touch the back of the foot and leg unit. *(Figure C)*
- The leg rest will lift for a few seconds and then automatically extend to fit your height. To make sure it adjusts correctly, place your foot fully back so your heel touches both the back and the bottom of the footbed.

Figure A

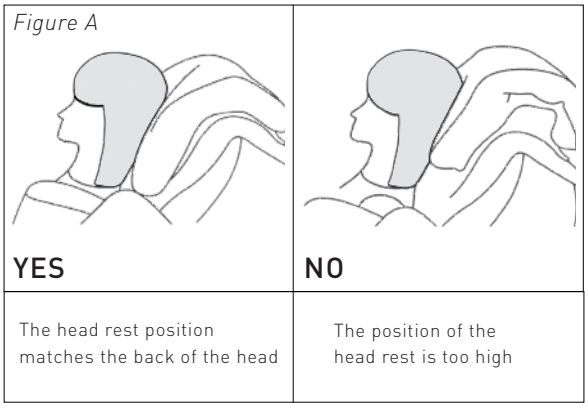


Figure B


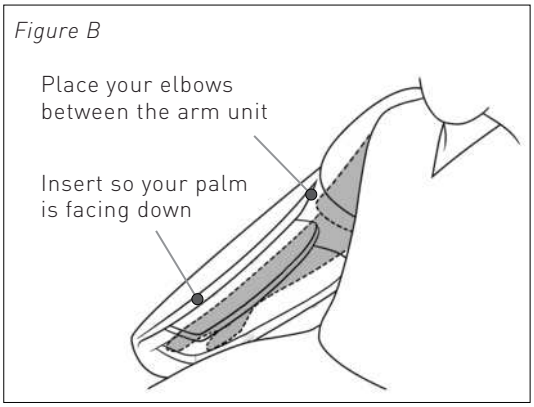


Figure C

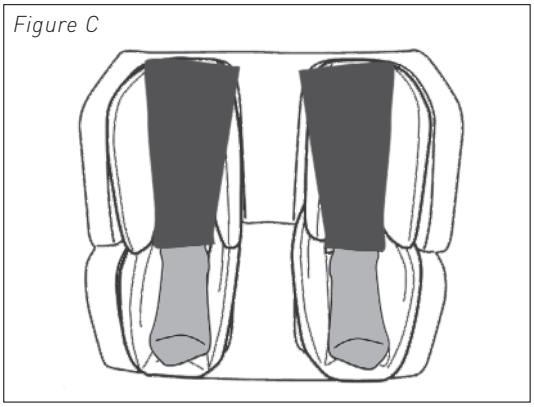

CAUTION

TO AVOID INJURY

- If you feel the massage is too strong, reduce the stimulation using the Massage Intensity button as described in this manual. If the massage remains too strong, stop operation immediately.
- Do not place hands or feet into the roller tracks, into the space between the seat and back, or into the moving parts.
- Do not sit wearing hard objects, such as hair accessories or jewelry.
- Keep massage sessions to no more than 30 minutes in length. Do not exceed 5 minutes for an area-specific massage for any part of the body.

Retrieving & Storing the Remote

- The remote should be stored on the remote stand when not using the chair.
- It is recommended to leave the remote on the remote stand during a massage. If choosing to retrieve the remote, gently pull it off of the remote stand.
- To return the remote to its storage position, gently place it on the stand and push until it snaps into the remote holder.

Powering On

- To turn the remote control power on, press the Power Button located on the face of the remote control.
- After powering on, the brand logo will briefly display before automatically transitioning to the Auto Programs screen.

Powering Off

- To turn the remote control power OFF, press and hold the Power Button located on the face of the remote screen.
- When the power is OFF, the remote screen is black/blank.
- If the chair has been in a reclined position it will automatically restore to the upright position when powered off.

NOTE: The Power Button on the remote is not the Main Power switch. See page 13 for Main Power.

Starting an Auto Massage

- Tap the "AUTO" button on the bottom of the screen. This will navigate you to the Auto Program menu, where you can select from a variety of pre-set massage programs designed for specific needs and preferences. Note, you can tap the forward/back arrows on the lower portion of the screen to select from additional program offerings. *(See Figure A)*
- Tap the desired massage button to activate that massage.
- Once the desired massage is active you will see the active massage summary/options display. *(See Figure B)*

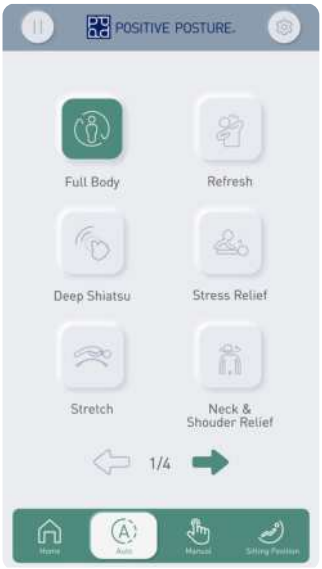


Figure A



Figure B

Shoulder Scan

To experience an auto massage as intended, it's important for the chair to accurately measure the height of your shoulders prior to the start of a programmed massage session. To ensure massage accuracy, please follow the steps below:

The massage chair will auto-detect your shoulder position by moving from your glutes to your shoulder top. During this time, your remote will display a screen indicating this process is in progress. *(See Figure A)*

Fine Tune Shoulder Position: Once auto-detect is complete, the remote will beep and a screen will appear that allows the roller position to be adjusted by tapping the Roller Up or Down arrows on the screen until the mechanism is located in the proper position at the tops of your shoulders. Once positioned correctly, wait 5 seconds for the scan position to register, or press OK. *(See Figure B)*

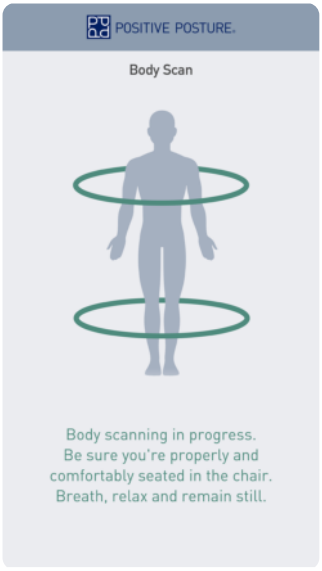


Figure A

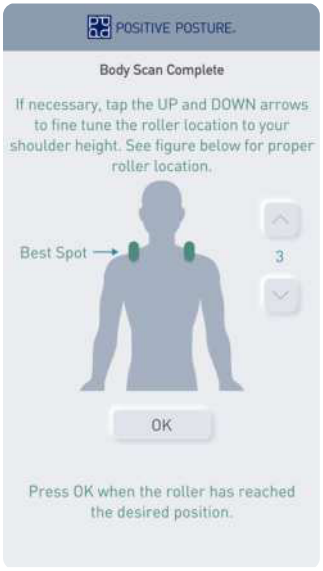


Figure B

Adjusting Auto Massage Settings

The Home screen appears after a massage program has been selected and allows for adjustments to enhance your massage experience.



1 Selected Program Indicator

The name of the currently active massage appears here.

2 Countdown Timer

A countdown timer displays the remaining time for your massage.

3 4D Intensity Adjustment

Tap the "4D Intensity" button, then tap the plus or minus button, to adjust the intensity, or the degree of pressure felt. The intensity can be adjusted between six levels, 1 being low and 6 being high. The default setting is 2.

4 Air Massage Intensity Adjustment

Tap the "Air Intensity" button, then tap the plus or minus button, to adjust the air intensity, or the degree of pressure felt from the air massage. The intensity can be adjusted between three levels, 1 being low and 3 being high. The default setting is 2.
Note: Air intensity adjustment is not available during Stretch programs.

5 CoreWave Intensity Adjustment


Tap the "CoreWave Intensity Adjustment" button to adjust the intensity, or the frequency of the CoreWave movement.
The intensity can be adjusted between three levels, 1 being low and 3 being high. The default setting is 2.

6 Time +/-

To change the length of your massage, tap the -5 or +5 buttons located on either side of the countdown timer. Each tap adjusts the session time by 5 minutes. You can set the timer for a minimum of 5 or up to a maximum of 30 minutes.

OPERATION & USE

Starting & Adjusting a Manual Massage

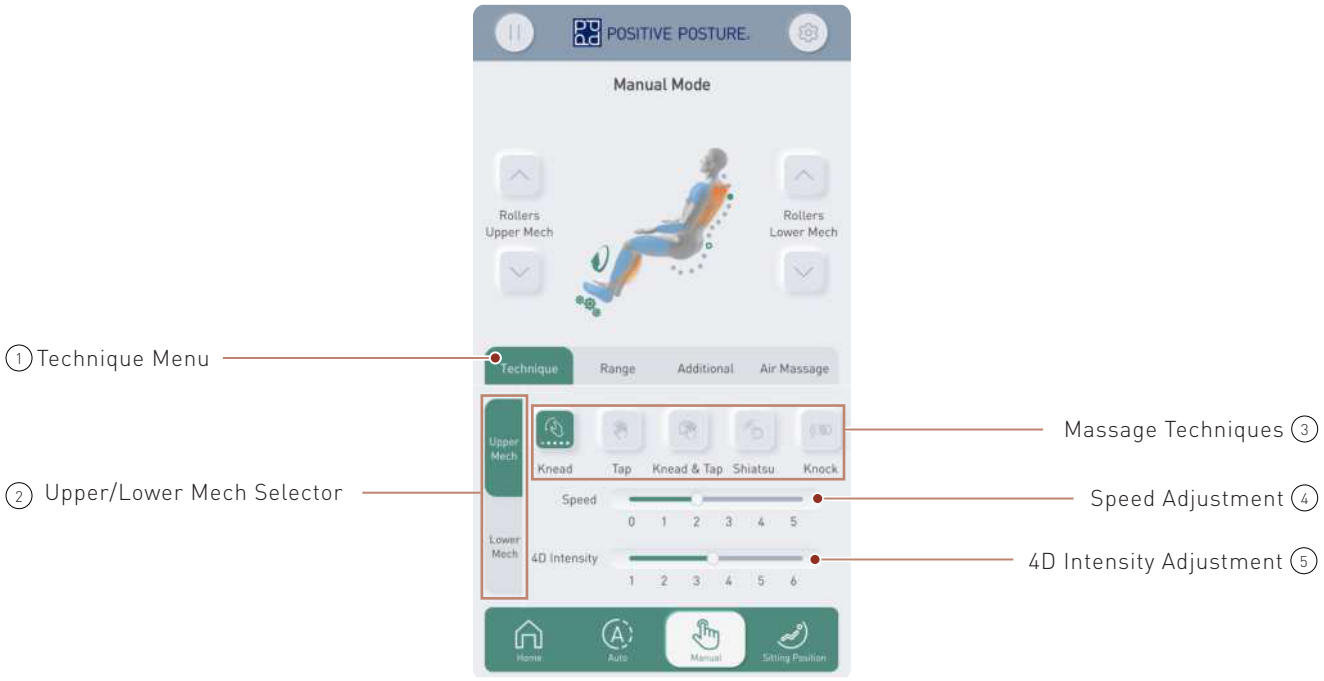

1 Technique Menu

Tap this menu button to view and select your preferred massage method.

2 Upper/Lower Mech Selector

Manual massages can be activated and customized for both the Upper and Lower mechanisms. Tap the "upper mech" or "lower mech" button, depending on which mechanism you wish to activate/adjust.

NOTE: 4D Intensity cannot be adjusted on the lower mechanisim.

3 Massage Techniques

Tap the desired massage method to start your session. For the upper mechanism, you can scroll left or right to see and select from additional options that do not fit on the main screen.

4 Speed Adjustment

Speed, or the pace of the massage mechanism movements, can be adjusted by moving the slider left or right to set your preferred setting. There are 5 speed levels, where 0 is off, 1 is the slowest, and 5 is the fastest.

5 4D Intensity Adjustment

4D Intensity, or the pressure of the massage mechanism, can be adjusted by moving the slider left or right to set your preferred setting. There are six levels, where 1 is the softest, and 5 is the firmest. The default setting is level 2.

NOTE: 4D Intensity cannot be adjusted on the lower mechanisim.

Starting & Adjusting a Manual Massage, continued...



1 Range Menu

Tap this menu button to adjust the roller's range of motion/position and width.

2 Rolling Range Adjustment

Rolling coverage, or the area of your back that is massaged, can be adjusted between Fixed Point, Partial, & Full Back by selecting the desired "Range" button. The default setting is Full Back. Feature available for all Manual methods.

NOTE: Mechanism will not move when the "Point" range is selected in the "Roll" massage mode.

3 Width Adjustment


Width, or the distance between the rollers of the massage mechanism, can be adjusted by moving the slider left or right between Narrow, Medium and Wide. The default setting is Wide. Feature available for the following Manual methods: Tap, Knock, Shiatsu, and Roll.

4 Using the Roller Position Arrows:

In Fixed Point, tap and hold the up/down Roller buttons to move the Upper or Lower massage mechanism to the desired location. The up arrow will move the rollers up the track toward the shoulders, while the down arrow will move the rollers down the track toward the glutes.

OPERATION & USE

Starting & Adjusting a Manual Massage, continued...

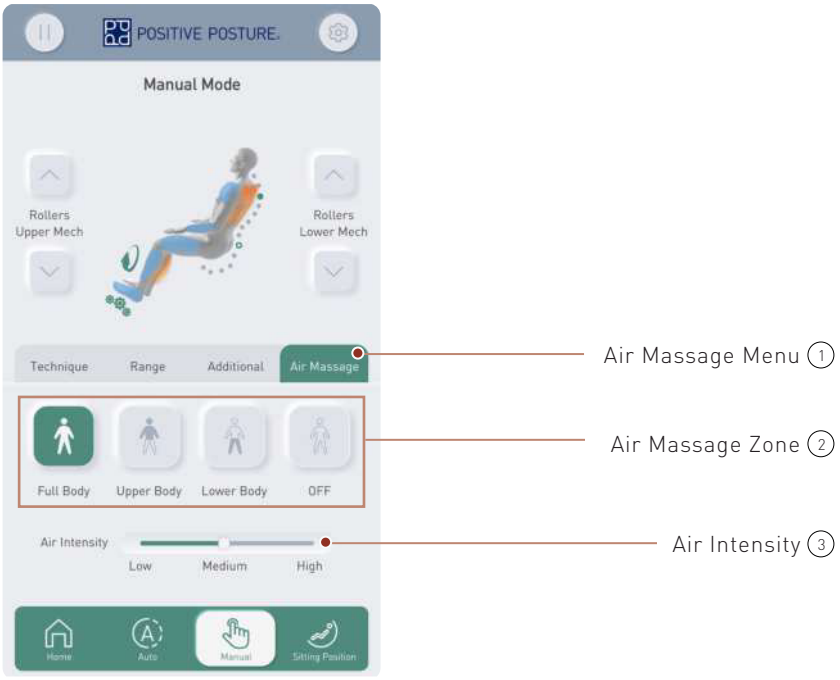

① **Air Massage Menu**

Tap this menu button to adjust the air massage's zone and intensity settings.

② **Air Massage Zone**

Tap the desired "ZONE" button to activate an Air Massage. Options include Full Body, Upper Body, Lower Body and Off.

③ **Air Massage Intensity Adjustment**

Air Intensity, or the pressure of the Air Massage , can be adjusted by moving the slider left or right to set your preferred setting. There are 3 levels, where Low is the lightest, and High is the firmest. The default setting is Medium.

Other Massage Settings

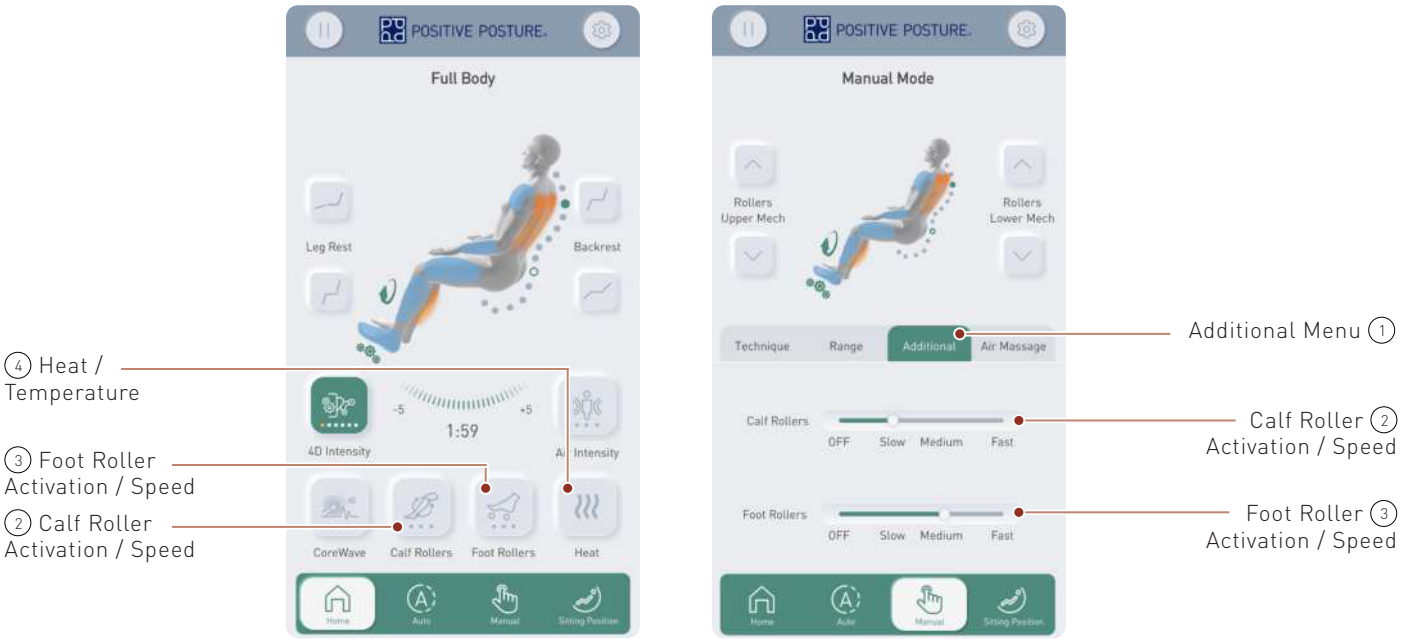


Figure A *Figure B*

① Additional Menu

Located under the Manual menu, tap this menu button to activate or adjust the speed of the calf & foot rollers. *See Figure B*

② Calf Roller Activation / Speed

When an Auto program starts, it automatically activates the calf rollers and TheraKnead™ Calf massage to run intermittently throughout the massage. The default speed is medium. To adjust the speed, or turn the Calf Rollers off, tap the "Calf Rollers" button located on the Home screen. *See Figure A*

Calf Rollers can also be activated as a manual massage by tapping the Additional Menu under the Manual Menu. The speed of the Calf Rollers can be adjusted by moving the slider left or right to set your preferred speed. There are 3 levels, Slow, Medium & Fast. *See Figure B*

NOTE: On the Home screen, highlighted Calf Rollers button indicates the feature is selected and active.

NOTE: Calf Rollers & TheraKnead™ Calf Massage will always activate together, they cannot be turned on and off seperately.

③ Foot Roller Use & Adjustment

Some of the Auto programs automatically activate the foot rollers to run intermittently throughout the massage. The default speed is medium. To adjust the speed, or turn the Foot Rollers off, tap the Foot Rollers button located on the Home screen. *See Figure A*

Foot Rollers can also be activated as a manual massage by tapping the Additional Menu from on the Manual menu. The speed of the Foot Rollers can be adjusted by moving the slider left or right to set your preferred speed. There are 3 levels, Slow, Medium & Fast. *See Figure B*

NOTE: A highlighted Foot Rollers button indicates the feature is selected and active.

④ Heat Activation / Temperature

On the Home page, tap the "Heat" button to turn the feature ON. A highlighted Heat button indicates the feature is active. Press once for low heat. Press again for medium heat, and a third time for high heat. Press once more to turn it off.

NOTE: All Auto programs automatically activate heat.

NOTE: Heat can be enjoyed on its own, or during any massage session.

OPERATION & USE

Chair Angle Adjustment

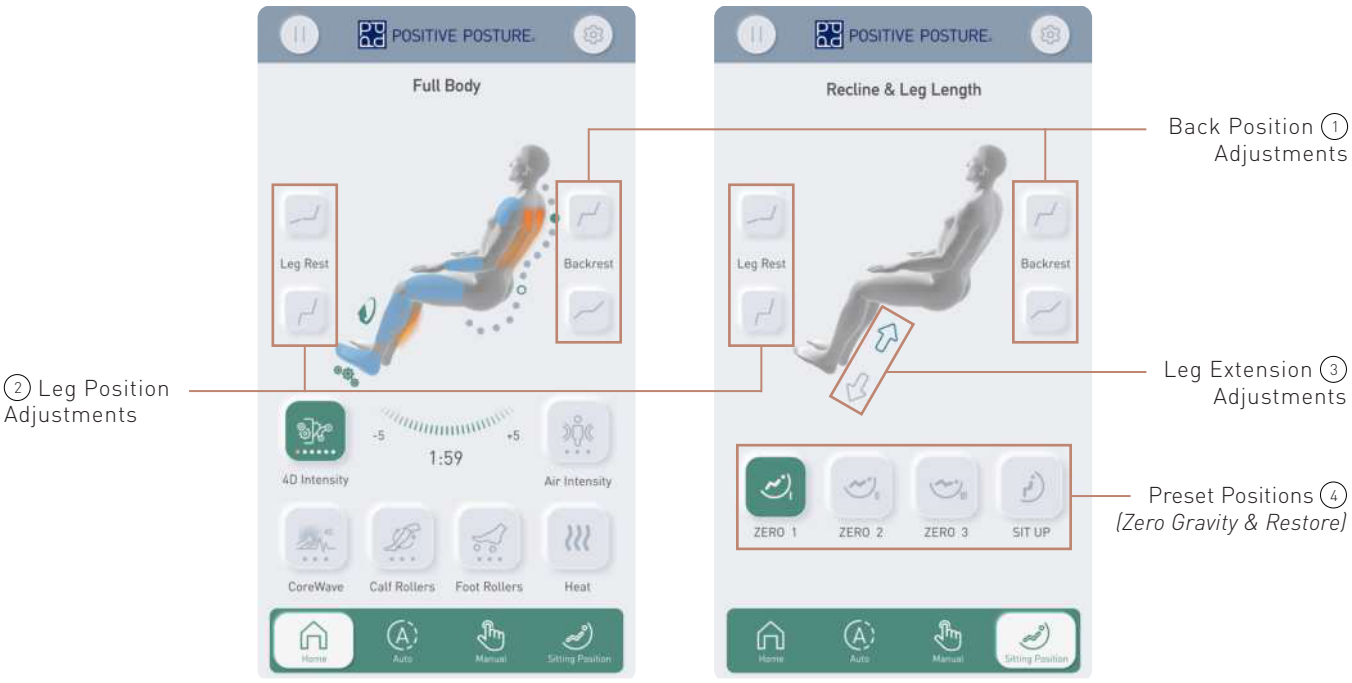


Figure A Home Screen *Figure B Sitting Position Menu*

During Auto massages, the chair will automatically recline to the suggested position for the massage. If a different reclining angle is desired you can adjust from the Home screen or the Sitting Position Menu as follows:

① **Back Position Adjustment**

To recline the angle of the chair's back, tap and hold the Backrest Down button until the desired position is reached. To incline the angle of the chair's back, tap and hold the Backrest Up button until the desired position is reached.

② **Leg Position Adjustment**

To recline the angle of the chair's leg, tap and hold the Leg rest Down button until the desired position is reached. To incline the angle of the chair's Leg rest, tap and hold the Leg rest Up button until the desired position is reached.

③ **Leg Extension Adjustment**

To extend the length of the leg unit, tap and hold the Down Arrow button located behind the leg of the human figure until the desired position is reached. To shorten the length of the leg unit, tap and hold the Up Arrow button located behind the leg of the human figure until the desired position is reached. *See Figure B*

④ **Pre-Set Positions (Zero Gravity & Restore)**

The chair is programmed with four pre-set positions that are designed for optimum comfort.

- Zero Gravity 1, 2, & 3: Tap the desired "Zero" button to quickly access the Zero Gravity position of your choosing. Note, the recline angle gets increasingly deeper from Zero 1 to Zero 2 to Zero 3.
- Restore: Tap "Restore" button to bring the chair to it's full and upright position.

Memory Settings

- To save time customizing settings exactly as you prefer them, the CoreWave allows you to save up to six favorite program and setting combinations (including recline position, 4D intensity, heat level, and air massage intensity). Once you ask the chair to memorize a Favorite, you can recall the selected message, and the CoreWave will perform the session just as you have customized it.



Setting a Favorite:

- Activate the desired massage and adjust the settings according to your preferences. When everything is set as desired, press the Auto Menu, then use the right Navigation arrow located at the bottom of the remote screen, to navigate to the Memory page. Press and hold the desired memory button for 3 seconds to set the Favorite session. When complete the icon will change from the standard Memory icon to the selected Auto Program Icon.

Running a Favorite:

- Press the Auto Menu, then use the right Navigation arrow located at the bottom of the remote screen, to navigate to the Memory page. Tap the desired Memory button to activate the session as previously saved.


Pausing a Massage

To pause a massage, tap the "Pause" button, on the top left side of the remote screen or press the physical pause button located on the lower right side of the remote. Upon doing this the massage mechanism will immediately cease all movement, the chair angle will remain in the exact position it was when the button was pressed. To resume the massage, tap either of the "Pause" buttons again. Upon doing this the massage will resume from where it was stopped.



Ending a Massage

To end a massage in the middle of an active session, press the "Power" button on the lower left of the remote. This will stop all massaging action, the rollers will move to their storage positions, and the chair will reset to it's fully upright position.

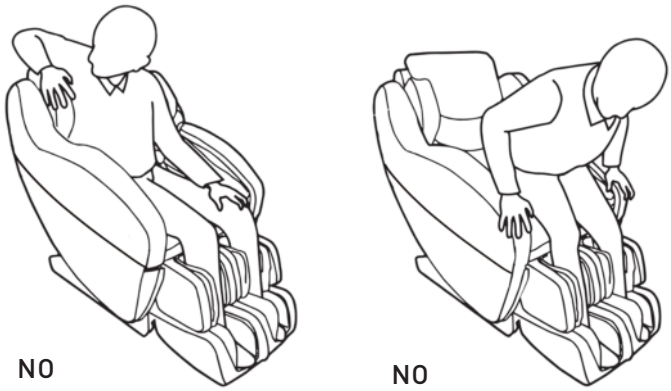

Remote Sound

The remote can be set up to either beep or remain silent with each button press. To change the Sound setting, tap the settings icon located on the top right side of the screen. This will navigate you to the Settings menu, where you can adjust Volume settings by moving the slider left or right to set your desired volume. 1 is quiet and 6 is the loudest. The default setting is OFF.

LCD Brightness

To adjust the brightness of the LCD display, tap the settings icon on the top right side of the screen. This will navigate you to the Settings menu, where you can adjust the screen brightness by moving the slider left or right to confirm your desired setting.

Speakers & Connecting to Bluetooth

- To use the speakers in the massage chair, connect a smart device to the chair's Bluetooth by doing the following:
- Ensure the Main Power is ON.
- From the settings menu of your smart device, select to connect to the Bluetooth connection titled POSITIVE POSTURE.
- Your Bluetooth device will indicate once the connection is established.
- Once a smart device is connected, use your smart device to select, play, pause audio and adjust the volume of sound.
- To connect a new device, the currently connected device must first be forgotten which can be done through the settings menu of your smart device.

Voice Control

The CoreWave is equipped with an optional Voice Control feature, which is set to OFF by default. To activate this feature, navigate to the Settings menu and select the ON button. Once enabled, you can wake the chair by saying "Hello, Positive Posture," and the chair will emit a soft chime to indicate it is ready for input. A complete list of recognized voice commands can be reviewed within the Settings menu by pressing the TIPS button next to Voice Control.

Simple View

Simple View optimizes the display for readability and ease of use by removing images and enlarging text and buttons. To turn Simple View ON, navigate to the Settings menu and press the icon next to Simple View. To exit Simple View and return to the standard display, touch the arrow button located on the top right of the screen.

Language

The remote control can be configured to display text in multiple languages. To change the display language, navigate to the Settings menu and select your desired language. All on-screen text and menus will update automatically to the selected language.

Screen Saver Display

The remote control features an adjustable screen saver function that automatically dims the display after a period of inactivity to encourage ultimate relaxation. The delay time can be customized between 30 seconds and 20 minutes. To adjust this setting, navigate to the Settings menu and move the slider left or right to set your deired Screen Saver Delay setting. Once the screen goes dark, simply tap anywhere on the display to reactivate it.

Accent Lighting (Exterior Armrests)

- Your chair features integrated accent lighting located on the exterior of the armrests. This subtle illumination enhances visibility and adds a touch of modern elegance.
- It is designed for seamless operation and requires no manual control. The default setting is directly linked to the chair's remote power system. When the chair's remote power is activated (i.e. the chair is on and responsive to the remote control), the accent lighting on the armrests will automatically illuminate. When the chair's remote power is turned off (i.e., the chair is off and unresponsive to the remote control), the accent lighting will automatically switch off.

Storage

- After every use, turn off the main power switch located on the back of the chair and unplug the power cord from the wall outlet.
- When not in use, the chair should always be stored in its full upright position (backrest fully inclined, footrest to its lowest position).
- The chair should be stored in a place of low humidity.
- Do not store the chair in direct sunlight or in places of high temperature. This may cause the color of the upholstery to fade.

Cleaning & Care



TO AVOID INJURY

- Before cleaning, always unplug the power cord from the wall outlet.
- Never clean the chair or remote with wet hands.



TO AVOID DAMAGE

- Before cleaning, always unplug the power cord from the wall outlet.
- Never clean the chair or remote with wet hands.

- Use a soft, dry cloth to wipe dust off of the remote control and the polyurethane surfaces located on the inside of the chair.

TROUBLESHOOTING

Troubleshooting Guide

If you encounter problems with the operation of your massage chair, please check the following points to see if they can be easily rectified. If the chair still malfunctions or does not work, please turn off the main power switch and remove the power plug from the wall outlet.

Contact Positive Posture for service. Do not dismantle or attempt to repair the chair yourself.
Doing so will void the warranty.

Problem: My chair doesn't work after powering on.
Possible Cause & Solution:

- Ensure the chair is securely plugged into the wall outlet.
- Ensure the main power switch on the back of the chair is switched to the ON position.
- Ensure the remote control and power cord are fully installed.
- Ensure the remote power is on.
- Ensure a massage session has been selected by touching any programmed or manual massage option on the remote.
- Other issues may include the wires or plug being damaged, the fuse having burned out, or an issue with the inner circuit. Call Positive Posture for service.

Problem: The power cord or plug is abnormally hot.
Possible Cause & Solution:

Stop using the chair immediately. Turn off the main power switch and remove the power plug from the wall outlet. Call Positive Posture for service.

Problem: My chair is very noisy.
Possible Cause & Solution:

There are many sounds made by the massage chair that are normal. This includes noise during CoreWave, tapping, the sound of the rollers rubbing against the seat covering while kneading, the sound of the rollers moving up and down the track, and the sound of the air pump and the air bags inflating and deflating. If you hear sounds other than this, or the sounds suddenly become different/louder, please call Positive Posture for service.

Problem: My chair makes a beeping sound at the start of a massage.
Possible Cause & Solution:

This sound is normal. There is no cause for concern.

Troubleshooting Guide, continued...

Problem: The rollers suddenly stopped in the middle of a session.
Possible Cause & Solution:

- The chair may have lost power. Turn off the main power switch, unplug from the wall, then reinsert the power cord into the wall outlet and turn on the main power switch.
- Ensure the massage session didn't complete its cycle by using the remote to select a new massage session.
- Sometimes when there is too much pressure on the rollers, they can stop moving temporarily. To test this, slightly lift your body or slightly raise the reclining seat back.

Problem: The body scan is not working properly.
Possible Cause & Solution:

- Ensure you are sitting still during the shoulder scan. If you wiggle or move around, you can trigger a false reading.
- Try sitting further back in the chair so that your head rests firmly on the head pillow. You can also adjust the shoulder position with the Roller buttons immediately after the scan is complete.

Problem: The backrest or leg rest will not go down entirely.
Possible Cause & Solution:

Ensure there are not any objects blocking its path.

Problem: The pillow is too high/too low.
Possible Cause & Solution:

Simply push up or pull down on the pillow to change its position.

Problem: The mechanism is not moving when the manual "Roll" massage mode is active.
Possible Cause & Solution:

- If the "Point" range is selected while the "Roll" massage mode is active, the mechanism will not move.
- Select a different range or massage mode to reactivate the mechanism.

Problem: The leg rest length adjustment is not working properly.
Possible Cause & Solution:

Ensure your feet are set so that your heel is touching the back of the footbed. If your feet are forward they will not trigger the sensor.

TROUBLESHOOTING

Troubleshooting Guide, continued...

Problem: My chair does not incline when I use the back up button.

Possible Cause & Solution:

- Press and hold the back up button, the chair should return to its upright position within this time.
- Turn off the main power switch and unplug the chair from the wall. Then re-insert the power cord and turn on the main power switch. The chair should return to its upright position.
- If the chair does not, turn the power on with the remote button, then turn the remote power off again. The chair should return to its upright position.
- If the problem persists, please contact Positive Posture service for further assistance.

Q: Can I put the chair on a wood floor?

A: Yes. However, because the chair could damage the floor surface, and because it will be louder, we recommend laying a carpet or something similar under the chair.

Q: I spilled water on the chair. What should I do?

A: There is danger of electric shock. Stop using the chair immediately and do the following:

- Turn off the main power switch on the back left of the chair.
- Unplug the power from the wall outlet.
- Contact Positive Posture for assistance.

Q: How do I make the massage stronger?

A: Try any of the following things:

- Adjust intensity to select maximum pressure.
- Sit back further in the chair so that your weight is against the reclining seat back.
- Lower the reclining seat back.

Q: How do I make the massage softer?

A: Try any of the following things:

- Adjust intensity to select minimum pressure.
- Raise the reclining seat back.

Warranty

We're confident that your new massage chair will provide you with many years of relaxation and trouble-free performance. Your Positive Posture massage chair includes a 3-year limited warranty. Some restrictions apply. See www.PositivePosture.com for full warranty details. We recommend that you print the limited warranty information at the time of purchase and keep it with your Owner's Manual.

Service & Repairs

Before calling for repairs or service, please consult the Troubleshooting Guide. If the problem persists, contact Positive Posture for assistance. For fastest service, please have the following information ready:

- Your name, email, address & telephone number
- Model Number (see the specification field on the front cover of this Owner's Manual or review the sticker affixed to the back base of your chair)
- Serial Number (see the sticker affixed to the back of your chair)
- Copy of original receipt and purchase date
- Malfunction or issue with as much detail as possible

SPECIFICATIONS

Model Name	CoreWave
Model Number	PP6218
Power Supply	AC110-220V | 50/60Hz
Power Consumption	240W
Rated Time	20 minutes
Preset Sessions	18 sessions
Manual Sessions	6 session on the upper mech, and 5 sessions on the lower mech. Foot Rollers, Calf Roller, and Full Body Air
Air Pressure	Approximately 38 kPa
Program Timer	Approximately 10-20 minutes
Footrest Adjustment Range	Approximately 5.5" (*14 cm*)
Reclining Angle	Approximately 128° to 156° against flat floor
Chair Dimensions	Upright: 65" *(165cm)* Lx 31" *(79cm)* W x 45.25 *(115cm*) H Reclined: 75.6" *(192cm)* Lx 31" *(79cm)* W x 33.5" *(85cm)* H
Seat Width	20.7" *(52.5cm)*
Seat Depth	19" *(48cm)*
Chair Weight	224 lbs *(102 kg)*
Supplied Items	Power cord, Remote Stand, Remote Control, Leg Rest
Shipping Dimensions	55" *(1410mm)* L x 31.9" *(810mm)* W x 47.2" *(1200mm)* H
Shipping Weight	257 lbs *(117 kg)*
Recommended User Weight	265 lbs *(120 kg)* (maximum)
Recommended User Height	5' 0" *(153 cm)* to 6' 4" *(193 cm)*

COREWAVE

by POSITIVE POSTURE®

Your massage chair was meticulously designed to conform to your specific comfort and relaxation needs. If you have questions about the operation of your massage chair, contact information for a Positive Posture chair expert can be found at www.PositivePosture.com.

Positive Posture. A division of Furniture For Life.

© 2026 Positive Posture

Boulder, Colorado | USA

www.PositivePosture.com